RESPIRATORY EQUIPMENT RENTALS

24 HOUR SERVICE

www.valuemed.co 888-998-8633 ValueMed.Sales@gmail.

ValueMed

YOUR JUST IN TIME RENTAL SUPPLIER

















CALL FOR IMMEDIATE DELIVERY


888-998-8633 Valuemed.co

INDEX

BIPAP
Philips V601
Respironics Vision
Philips V-305
System One5
VENTILATORS
Puritan Bennett 9802
Puritan Bennett 840
Maquet Servo I
Philips Trilogy 202
LTV 1200
MRI LTV 1200
WIN LI V 1200
HUMIDIFICATION
Fisher Paykel MR 850 Heater4
HIGH FLOW OXYGEN THERAPY
Airvo 26
Airvo 2
Airvo 2
Airvo 26
Airvo 2
Airvo 2
Airvo 2
Airvo 2
Airvo 2 6 Optiflow / Maxventuri 6 Vapotherm 7 High Flow Blender/Heater 7 THERAPY EQUIPMENT The Vest 8 Metaneb 8
Airvo 2
Airvo 2 6 Optiflow / Maxventuri 6 Vapotherm 7 High Flow Blender/Heater 7 THERAPY EQUIPMENT The Vest 8 Metaneb 8
Airvo 2



DEADIDATABY EQUIDATENT	
RESPIRATORY EQUIPMENT	
SALES	
Oxygen Analyzers	12
Percussor	
Respirometer	12
Space Saving IV Pole	
Suction Regulators	
Hoses & DISS Adapters	
1	
BLENDERS & ACCESORIES	
Max Venturi	13
Brackets	13
Flowmeters	13
INFANT CARE	
Infant Flow SiPAP	14
Infant Flow NCPAP	14
VALUEMED PROFESSIONAL	
SERVICES	

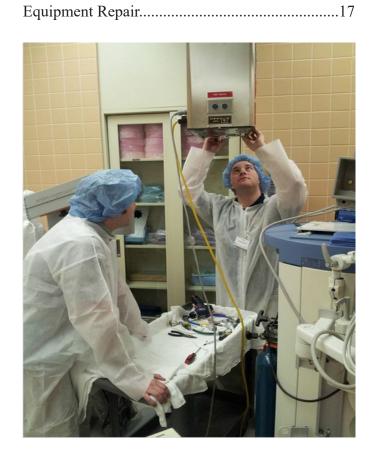



24 HOUR SERVICE
YOUR JUST IN TIME RENTAL SUPPLIER

PHILIPS V-60			
FOR USE WITH PART NUMBER		DESCRIPTION	
	RENTAL-RM-801	Fisher Paykel Heater and Cable/Probe	
	RENTAL-RM-505	Transport Oxygen Manifold	
V-60	V-60-6F	5 Pack of Inlet Filters	
	RB-002	Support Arm Rail Bracket	
	116-56	Support Arm	

The Philips Respironics V60 Ventilator is an assist ventilator that is intended to augment patient breathing. It is intended for spontaneously breathing individuals who require mechanical ventilation.



V60 AVAILABLE WITH HIGH FLOW THERAPY CAPABILITY

-Allows clinician to switch between NIV and HFT as patient condition changes

HEATER & CABLES **AVAILABLE**



RESPIRONICS BIPAP VISION		
FOR USE WITH PART NUMBER		DESCRIPTION
	RENTAL-RM-801	Fisher Paykel Heater and Cable/Probe
BiPAP VISION	VIS-6F	Replacement Inlet Filter 6/pack
DITAL VISION	RB-002	Support Arm Rail Bracket
	116-56	Support Arm

The BiPAP Vision System offers a form of sub-acute therapy for patients and may eliminate the need for intubation in appropriate candidates. The integrated display screen uses easy-to-read real time graphics in waveforms or bar scale format to provide enhanced monitoring for caregivers. Auto-Trak Sensitivity ensures optimum sensitivity despite changing breathing patterns or circuit leaks.





VENTILATORS



PURITAN BENNETT 980 VENTILATOR			
FOR USE WITH PART NUMBER DESCRIPTION			
DD 000	RENTAL-RM-801	Fisher Paykel Heater With Cable/Probe	
PB 980	10043551	Expiratory Filter	

The Puritan Bennett 980 ventilator system helps patients breathe more naturally through innovative breath delivery technology. The Puritan Bennett 980 has an intuitive and highly configurable touch screen display that gives clinicians the options to see as little or as much information as needed including five graph and loop layout options. The 980 is equipped with an expiratory filtration system. The ventilator assurance feature will, in the event of certain system failures, will continue to deliver ventilatory support as close to the preset settings as feasible.

HEATER & CABLES **AVAILABLE**





PURITAN BENNETT 840 VENTILATOR		
FOR USE WITH PART NUMBER DESCRIPTION		
	RENTAL-RM-801	Fisher Paykel Heater With Cable/Probe
PB 840	10043551	Expiratory Filter
	Max-16	Replacement Oxygen Sensor

The Puritan Bennett 840 Ventilator System is A critical care ventilator. It is highly responsive and offers superior comfort, delivering sensitive, precise breaths to critically ill neonatal through adult patients.

All units supplied with air compressor

NEO AND ADULT MODES AVAILABLE



24 HOUR SERVICE

YOUR JUST IN TIME RENTAL SUPPLIER


VENTILATORS

MAQUET SERVO-I VENTILATOR			
FOR USE WITH PART NUMBER DESCRIPTION			
SERVO-I	RENTAL-RM-801	Fisher Paykel Heater With Cable/Ptobe	
SERVO-1	64-87-180	Servo I - Replacement battery	

The Servo – I ventilator is a flexible ventilation platform and is available for both infants and adults. These configurations for the Servo-I are equipped with a number of ventilation modes that are suitable for each patient category. Ventilation is accomplished via either Invasive or Non-Invasive Intubation. The Non-Intubation mode also features leakage compensation, the Servo-I also has a transport mode.



HEATER & CABLES **AVAILABLE**

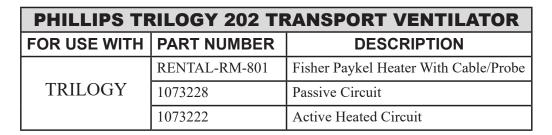

With one simple setting change, the Trilogy 202 supports either active or passive exhalation breathing circuits to accommodate changes in circuit preference. Compact design with long-life internal battery for intra-hospital transport. Supports active and passive circuits for invasive and noninvasive ventilation

BATTERY FOR TRANSPORT







VENTILATORS



LTV 1200 & ** MRI LTV 1200 TRANSPORT VENTILATOR**			
FOR USE WITH PART NUMBER DESCRIPTION			
	RENTAL-RM-801	Fisher Paykel Heater With Cable/Probe	
LTV 1200	29657-001	Disposable Adult Circuit Without PEEP	
L1 V 1200	29665-001	15 foot MRI Circuit Without PEEP	
	LTV-DT-PM	Pole Mount Bracket for LTV Ventilator	

The small, lightweight Vyaire LTV 1200 Ventilator provides sophisticated critical care ventilation in austere environments. The LTV Ventilator supports both invasive and noninvasive applications. It combines portability, impact-resistant construction, easy-to-read displays, and comprehensive monitoring. A single solution for any patient weighing 5 kg or more

SPRINT BATTERY PACK INCLUDED



MRI UNIT AVAILABLE



HEATER & CABLES AVAILABLE

FISHER PAYKEL MR850 HEATER AND CABLES		
FOR USE WITH PART NUMBER DESCRIPTION		DESCRIPTION
MR-850	RENTAL-RM-801	Fisher Paykel Heater With Cable/Probe
MIK-03U	RT202	Circuit & Chamber Kit
MR-850	BR-825-007	Pole Mount Bracket



One device supports multiple therapies for adult, pediatric and neonatal patients. Effective. One touch therapy selection combined with dual temperature feedback systems ensure ease of use. The MR850 uses flow sensor technology to automatically standby during periods of gas flow interruption. The MR-850 uses automatic algorithms for condensation control and humidity compensation.



24 HOUR SERVICE
YOUR JUST IN TIME RENTAL SUPPLIER

BIPAP

PHILLIPS V-30		
FOR USE WITH PART NUMBER DESCRIPTION		
V 30	RENTAL-RM-801	Fisher Paykel Heater With Cable/Probe

The Respironics V30 ventilator Auto BiPAP airway management system offers cost-effective treatment of obstructive sleep apnea patients. The V30 Auto features:

Nine flexible ventilation therapy modes
Auto-titrating BiPAP
Customizable ventilator alarms, a simple interface, and a convenient roll-stand

The V30 Auto is designed to improve operational efficiency while addressing the needs of your low-acuity, noninvasive ventilation patients.



RESPIRONICS SYSTEM ONE BIPAP AUTO		
FOR USE WITH PART NUMBER DESCRIPTION		
SYSTEM ONE	1063785	Disposable Water Chamber

Philips Respironics revolutionized sleep therapy by introducing bi-level positive airway pressure technology to treat obstructive sleep apnea. Emulating natural breathing with pressure relief at inhalation and exhalation, this proven technology does more than provide highly effective therapy – it delivers a level of comfort that can mean the difference between compliance and non-compliance for patients having difficulty adapting to traditional CPAP therapy.





24 HOUR SERVICE YOUR JUST IN TIME RENTAL SUPPLIER

HIGH FLOW OXYGEN THERAPY



FISHER PAYKEL AIRVO 2		
FOR USE WITH	PART NUMBER	DESCRIPTION
	RENTAL-RM-801	Fisher Paykel Heater With Cables
	900PT600	Red Disinfection Tube
AIRVO 2	900PT501	Tube and Chamber Kit
	OPT942	Small Nasal Cannula
	OPT944	Medium Nasal Cannula
	OPT946	Large Nasal Cannula

The AIRVO 2 system features a humidifier with integrated flow source that delivers high flows of air/oxygen mixtures to spontaneously breathing patients, through a variety of interfaces up to 60 L/min. The AIRVO 2 is an integrated system providing versatility, mobility and convenience.

IDEAL FOR ROOMS WITHOUT PIPED MEDICAL AIR



FISHER PAYKEL OPTIFLOW / MAXVENTURI			
FOR USE WITH PART NUMBER DESCRIP		DESCRIPTION	
	RT202	Circuit & Chamber Kit	
	OPT942	Small Nasal Cannula	
OPTIFLOW / MAXVENTURI	OPT944	Medium Nasal Cannula	
MAXVENTUR	OPT946	Large Nasal Cannula	
	MAX250-E	Replacement Oxygen Sensor	

The MaxVenturi has been specifically designed to work with high-flow Optiflow canulas. Using the Venturi principle, the MaxVenturi allows users to mix ambient air and oxygen without the need for piped medical-grade air. This is the perfect solution for adult high-flow oxygen therapy in almost any environment. Maxventuri delivers up to 60 L/Min of blended gas.

IDEAL FOR ROOMS WITHOUT PIPED MEDICAL AIR



24 HOUR SERVICE
YOUR JUST IN TIME RENTAL SUPPLIER

HIGH FLOW OXYGEN THERAPY

VAPOTHERM-PRECISION FLOW			
FOR USE WITH PART NUMBER DESCRIPTION		DESCRIPTION	
VAPOTHERM	PF-DPC-HIGH	High Flow Kit-Circuit & Cannula 5-40 LPM	
PRECISION	PF-DPC-LOW	Low Flow Kit-Circuit & Cannula 1-8 LPM	
FLOW	MAX250	Replacement Oxygen Sensor	

Vapotherm Precision Flow is designed for ease of use in mind. Control flow, oxygen, and temperature independently with the push of a single button. Vapotherm can be used on Neonates to Adults.

Flow up to 40 LPM



HIGH FLOW BLENDER AND FISHER PAYKEL HEATER			
FOR USE WITH	PART NUMBER	DESCRIPTION	
	RT202	Circuit & Chamber Kit	
HIGH FLOW BLENDER &	OPT942	Small Nasal Cannula	
HEATER	OPT944	Medium Nasal Cannula	
	OPT946	Large Nasal Cannula	

The High Flow Air/Oxygen Blender / Optiflow system uses a blender, 70 and 15 LPM flowmeter to deliver precise, humidified, blended gasses to patients up to 70 LPM. Oxygen concentrations are adjustable form 21% to 100%. When used in conjunction with the Optiflow nasal canulas, this system offers the ability to comfortably deliver a wide range of oxygen therapy.






24 HOUR SERVICE
YOUR JUST IN TIME RENTAL SUPPLIER

THERAPY EQUIPMENT

HILL-ROM THE VEST			
FOR USE WITH PART NUMBER		DESCRIPTION	
	P300629000	X-Small Wrap - 19 inches-23 inches	
THE VEST	P300630000	Small Wrap - 23 inches-33 inches	
UNIT	P300631000	Medium Wrap - 33 inches-43 inches	
UNII	P300632000	Large Wrap - 43 inches-53 inches	
	P300633000	X Large Wrap - 53 inches-67 inches	
	P300634000	XX Large Wrap - 67 inches-75 inches	

The Vest Airway Clearance System is designed to assist patients in the mobilization of retained secretions. If not removed, retained secretions may contribute to increased rates of respiratory infection, hospitalization, and reduced lung function.

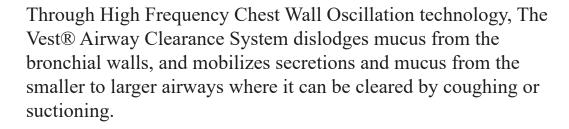

HILL-ROM METANEB

The MetaNeb System is indicated for mobilization of secretions, lung expansion therapy, the treatment and prevention of pulmonary atelectasis, and also has the ability to provide supplemental oxygen when used with compressed oxygen.

Metaneb introduces continuous high frequency oscillation (CHFO) therapy immediately after lungs have been expanded through con-

tinuous positive expiratory pressure (CPEP) The Metaneb also integrates aerosol delivery of medications for maximum efficiency.







THERAPY EQUIPMENT

RESPIRONICS COUGH ASSIST

Cough Assist safely removes broncho-pulmonary secretions by stimulating individual to cough. This is a desired approach to respiratory therapy as it's less invasive over suctioning the bronchial secretions. Respiratory Therapy is vital for individuals with respiratory conditions that make it difficult to cough and therefore mucus builds up which increases risk of respiratory infections. This model offers both automatic and manual modes of operation. Therapy involves slow inhale and rapid exhale of airflow to create positive to negative pressure. This loosens the mucus to stimulate a cough consisting of five coughing cycles.



PHILIPS RESPIRONICS COUGH ASSIST T-70

CoughAssist T70 is a mechanical insufflation-exsufflation device designed to noninvasively clear secretions from the lungs by simulating a natural cough. Similar to a normal deep breath, CoughAssist gradually applies positive air pressure (insufflation) to obtain a large volume of air within the lungs. The device then quickly reverses the flow of air by shifting to negative air pressure

The resulting high expiratory flow helps mobilize secretions out of the airway as a deep, natural cough would do.





SUCTION EQUIPMENT & SERVICES

LET US HELP WITH YOUR EMERGENCY AND PLANNED SHUTDOWNS * PLANNING * SETUP * MONITORING * REMOVAL * COMPLETE SUPPORT OF SHUTDOWN * FULL SERVICE SHUTDOWN



MOBLVAC 3			
FOR USE WITH	PART NUMBER	DESCRIPTION	
Suction Pump	RW-067	Patient Ready Canister and Tubing Kit	

2 MODELS AVAILABLE

Continuous / Intermittant model- up to 380mmHg & 36 lpm flow (constant mode) 8 lpm (intermittent mode) IDEAL FOR FLOOR UNITS AND ICU SUPPORT

Continuous only model- up to 550 mmHg & 38 lpm flow IDEAL FOR SURGERY



AEROS INSTAVAC II

The Instavac II is a high vacuum & high flow suction unit. It produces constant suction with vacuum & flow surpassing standard wall suction. Constant suction providing 550 mmHg and over 38 lpm flow.

SSCOR DUET BATTERY OPERATED SUCTION PUMP



- An AC suction unit, with battery back up, the Duet provides instant, powerful suction whether or not AC power is available.
- When fully charged, the sealed lead acid battery provides 45 minutes of suctioning power.
- •Provides suction up to 525mmHg ***CRASH CARTS**



24 HOUR SERVICE
YOUR JUST IN TIME RENTAL SUPPLIER

SHUTDOWN EQUIPMENT & SERVICES

LET US HELP WITH YOUR EMERGENCY AND PLANNED SHUTDOWNS

* PLANNING * SETUP * MONITORING * REMOVAL * COMPLETE SUPPORT OF SHUTDOWN * FULL SERVICE SHUTDOWN

SHUTDOWN REGULATORS AND HOSES

- AIR - OXYGEN - NITROUS OXIDE - NITROGEN

50 PSI PRESET REGULATORS











CYLINDER CARTS FOR GAS SHUTDOWNS

Over 300 cylinder carts available in various tank sizes









AIR COMPRESSORS



Standard Compressor

Timeter 3500 42 LPM Flow 50 PSI Heavy Duty Compressor

Timeter 2000 56LPM Flow 50 PSI

Ideal for use with Ventilators





24 HOUR SERVICE
YOUR JUST IN TIME RENTAL SUPPLIER

RESPIRATORY EQUIPMENT & SALES

OXYGEN MONITORS & ANALYZERS











Cufflator

Repair

Services










W/Sensor IRC-10LX-O2 5L PerfectO2 W/Sensor **IRC-5PW**



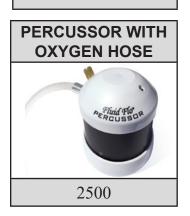









BLENDERS AND ACCESSORIES





MAX-VENTURI WITH SENSOR

MAX-250E Sensor



CUSTOM CONFIGURATIONS Elbowed Flowmeters, Silencers, Brackets, Hoses

Blender Accessories

10 FT Oxygen Hose, DISS Male x DISS Female BM-108-007 10 FT Air Hose, DISS Female x DISS Female BM-208-007 Elbowed 0-15 LPM Oxygen Flowmeter BM-300-007 Elbowed 0-3 LPM Oxygen Flowmeter BM-400-007 Elbowed 0-1 LPM Oxygen Flowmeter BM-500-007 Elbowed 0-70 LPM Oxygen Flowmeter BM-570-007 Elbowed 15 LPM Oxygen Flowmeter Wye BM-600-007

Air Inlet Filter/Water Trap AFA11629 Stand with Wheels 1059-007

Brackets

Pole Mount Bracket BR-800-007 Female Dovetail Bracket BR-825-007 Wall Mount Bracket BR-850-007

Hill Rom / Fairfield Blender Bracket BR-901-007 Suction Slide Bracket to Blender BR-903-007 Vertical Rail Bracket, GE Incubator BR-904-007

MAXVENTURI MUFFLER

R211P38



10 FOOT AIR / OXYGEN BLENDER HOSES


BM-208-007




BM-108-007

DAMPENS NOISE AND MINIMIZES PATIENT STIMULATION

HIGH FLOW SILENCER FITS 15, 22MM HOSE & FLOWMETERS 219P50

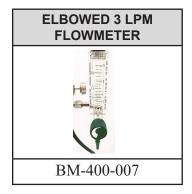



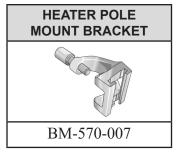






24 HOUR SERVICE YOUR JUST IN TIME RENTAL SUPPLIER

FLOWMETERS

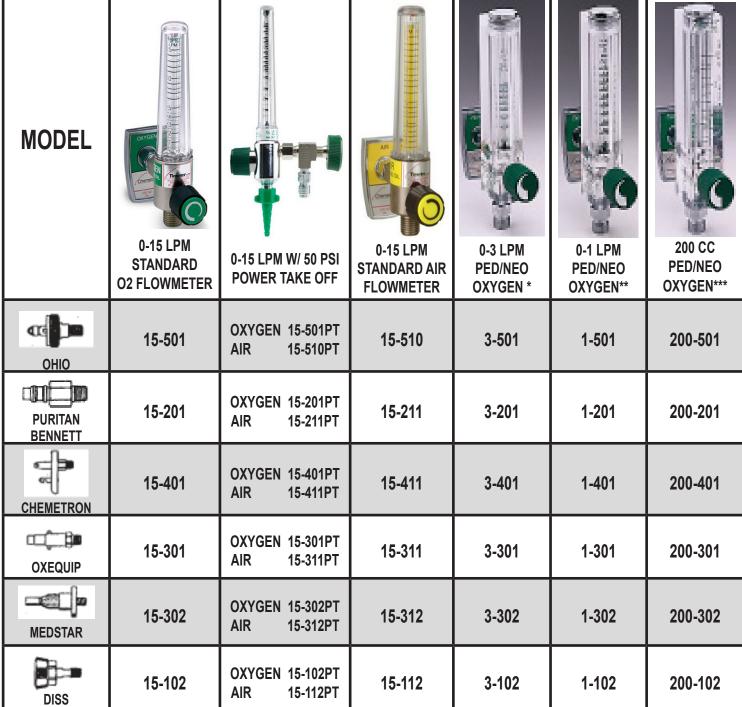



22-501-500	OXYGEN	\$144.00
22-511-510	AIR	\$149.00
22-531-530	VACUUM	\$169.00



CHEME	TRON WYE ASSE	MBLY
	0	
22-401-400	OXYGEN	\$144.00
22-411-410	AIR	\$149.00
22-431-430	VACUUM	\$169.00

MEDSTAR WYE ASSEMBLY		
		VOEW CO
22-302-300	OXYGEN	\$144.00
22-312-310	AIR	\$149.00
22 332 330	VACIIIIM	\$160.00



SUCTION EQUIPMENT

SUCTION REGULATORS













Ohio Int/Cont

Ohio Push to set

Ohio Continuous

Gentec

Precision Medical

SUCTION REGULATOR QUICK CONNECT BACK / WALL CONNECTOR OPTIONS

DISS FEMALE













PURITAN



SUCTION REGULATOR BOTTOM CONNECTOR OPTIONS



1/8" NPT DISS Male, Vacuum 25-132-02



Tubing Nipple Vacuum Regulator 25-999







130-1000 295-1000

VENTILATOR / BIPAP SUPPORT ARM AND BRACKET

116-56 **RB-002** *FLAT RATE SUCTION **REGULATOR REPAIR*** \$89.00

ALL MODELS * Chemtron * Ohmeda * Precision Medical * * Boehringer * Gentec *



24 HOUR SERVICE YOUR JUST IN TIME RENTAL SUPPLIER

HIGH PRESSURE HOSES & DISS ADAPTERS

HOSE QUICK CONNECT ADAPTERS



DISS PURITAN	PART #	GAS
BENNETT	201-55	O2
	211-55	AIR
	221-55	N20
	231-55	VAC
	251-55	EVAC
	281-55	CO2

DISS MEDSTAR	PART#	GAS
	302-55	O2
Gr	312-55	AIR
	322-55	N20
	332-55	VAC
- Call	354-55	EVAC
	355-55	CO2

DISS OHIO	PART#	GAS
	501-55	O2
	511-55	AIR
	521-55	N20
	531-55	VAC
	551-55	EVAC
	581-55	CO2

HIGH PRESSURE HOSES-BY THE FOOT / CUSTOM CONFIGURATIONS



- *Available in custom lengths and configurations
- *Available in a variety of gasses
- *Sold by The Foot *Built to your needs*
- *200 PSI Max

STANDARD 10 FOOT DISS FEMALE X DISS FEMALE HOSES

10 FT CONDUCTIVE OR / ANESTHESIA			
#100-10-22-100	Green	Oxygen	
#110-10-23-110	Yellow	Air	
#120-10-24-120		Nitrous Oxide	
130-10-25-130	White	Vacuum	
#140-10-26-140	Black	Nitrogen	
150-10-27-150	Lavender	EVAC	
180-10-28-180	Gray	Carbon Dioxide	

10 FT NON CONDUCTIVE RESPIRATORY		
100-10-20-100	Green	Oxygen
110-10-21-110	Yellow	Air











24 HOUR SERVICE
YOUR JUST IN TIME RENTAL SUPPLIER

PROFESSIONAL SERVICES

MED-GAS VERIFICATION SERVICES

LET US HELP WITH YOUR THIRD PARTY MEDICAL GAS VERIFICATIONS

ZONE VALVE BOXES



Dew Point Readings
Hydrocarbon Sampling
Purity Sampling
Trace Gas Sampling

ALARM PANELS



GAS OUTLET STATIONS
OUTLETS AND PARTS AVAILABLE







Regulators

ON SITE PURITY TESTING & RESULTS



ANNUAL TESTING
AFTER RENOVATION TESTING
INSTANT REPORTS TO
ENSURE SAFETY

EQUIPMENT REPAIR SERVICES

Flowmeters



Percussors

Wyes



Suction Regulators-Flat Rate Repair

Respirometers



Blenders





24 HOUR SERVICE
YOUR JUST IN TIME RENTAL SUPPLIER

ValueMed

YOUR JUST IN TIME RENTAL SUPPLIER

MEDICAL GAS SHUTDOWN SERVICES

LONG TERM & ON DEMAND RENTALS

MEDICAL GAS TESTING AND CERTIFICATIONS

SALESREPAIRS**





24 HOUR SERVICE